



HOCT-1/1F

All-in-One Optical Coherence Tomography with Fundus and Angiography

Specification

Type	Spectral Domain OCT, Fundus Digital Photography
Resolution in Tissue	20um(Lateral), 6-7um(z-axis) at index 1.36
Scan Speed	Max. 68,000 A-scans/sec.
Scan Range	X: 6-12mm, Y: 6-9mm, Z: 2.34mm
Acquisition Time of 3D Image	1.4 sec.(Normal Mode, A512 x B96)
Minimum Pupil Diameter	2.5mm
Light Source	840nm
Optical Power at Cornea	≤650mW
OCT	<p>Macular: Macular Line, Macular Cross, Macular Radial, Macular 3D, Macular Raster</p> <p>Disc: Disc Circle, Disc Radial, Disc 3D, Disc Raster</p> <p>Angiography Range X: 3-12mm, Y: 3-9mm(Optional)</p> <p>Angiography Map Superficial, Deep, Outer, Choriocapillary, Retina, Custom, Enface, Thickness, Depth coded Map</p> <p>Angiography Analysis FAZ, Vessel density Tools(Optional)</p> <p>Camera Built-in 12 Megapixel Color</p> <p>Angle of View 45°</p> <p>Minimum Pupil Diameter 4.0mm(Normal Mode), 3.3mm(Small Pupil Model)</p>
Fundus	<p>Flash Light White Light, 10 Levels</p> <p>Resolution Center: 60 Line pair/mm Middle: 40 Line pair/mm Periphery: 15 line pair/mm</p> <p>Working Distance 33mm</p> <p>Display 12.1 inch, 1280x800 pixel, Touch Panel Color LCD</p> <p>Dioptric Compensation -33D~+33D Total -13D~+13D with No Compensation Lens +7D~+33D with Plus Compensation Lens -33D~-7D with Minus Compensation Lens</p>
Common	<p>Internal Fixation Target LCD(internal), White LED(external)</p> <p>Horizontal Movement 70mm(back and forth), 100mm(left and right)</p> <p>Vertical Movement 30mm</p> <p>Chinrest Movement 62mm(up and down), Motorized</p> <p>Auto Tracking 30mm(up and down) 10mm(right and left) 10mm(back and forth)</p> <p>Network DICOM File support(need to be customized) Built in Computer</p>
Power Supply	AC100-240V, 50/60Hz, 16-07A
Dimensions / Mass	330(W)x420(D)x529(H)mm / 30kg
Accessories(Optional)	
Anterior Segment Module	<p>Scan Pattern ACA Line, Anterior Radial</p> <p>Scan Range 6-9mm(width), 2.3mm(depth)</p> <p>Software Analysis Corneal Layers, Thickness Map, Thickness & Angle</p>
Wide Anterior Segment Module	<p>Scan Pattern ACA Line, Anterior Radial</p> <p>Scan Range 16mm(width), 2.3mm(depth)</p>
Web-Viewer	Web-based, Multi users can be accessible. Progression Analysis, Comparison Analysis, 3D Analysis

Designs and details above can be changed without prior notice for the purposes of improvement.

Huvitz

TECNOIMAGEN
www.tecnoimagen.com.ar | 0810 333 8273

Innovative
Ophthalmology
Solutions

All-in-One Optical Coherence Tomography with Fundus and Angiography

HOCT-1/1F



Huvitz *redefine. Re^create*




3D OCT FUNDUS ANGIOGRAPHY

Huvitz All-in-One OCT

Combination of Fundus Camera and Angiography in 3D OCT

All-in-One HOCT is smart.
3D OCT Angiography and Fundus Camera, Totally integrated system combined with PC.
All-in-One HOCT is easy to use.
One button, it creates a High-speed Scan and a High-quality Image, giving greater insight for the ophthalmology clinic.
One space, allows the acquisition of various types of analysis and diagnosis results.
Easy to use and the results are outstanding and easily follow.
Huvitz All-in-One HOCT will be the icon for leading a new era of Optical Coherence Tomography(OCT).

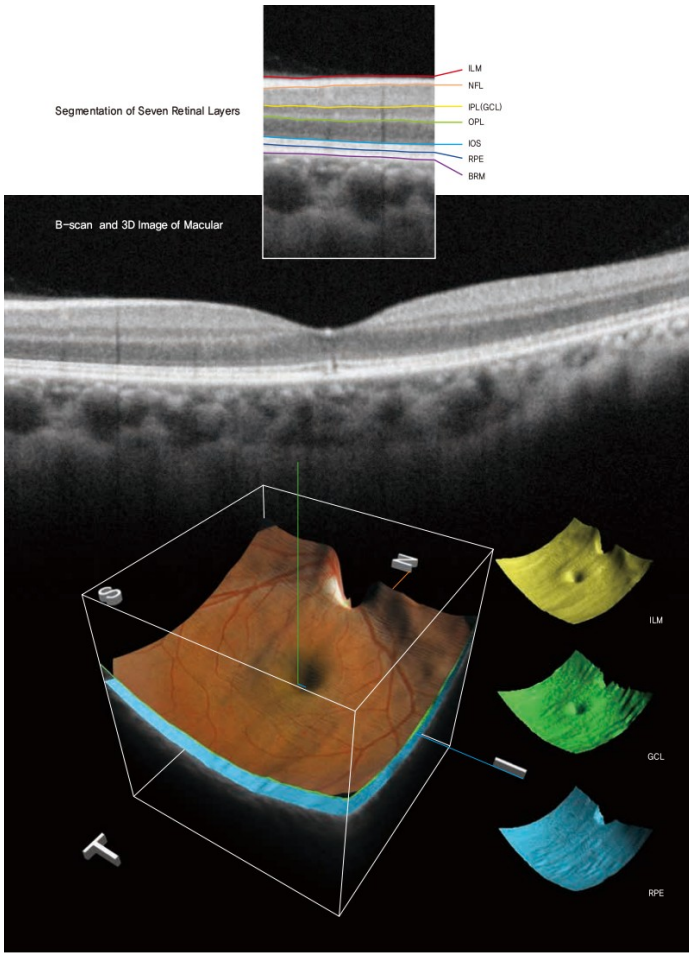

High-Speed & High-Quality

Incredible speed of 68,000 A-scan/sec. :
More Realistic and Clearer image in high resolution

Provides High-speed Scan, High-quality Image by using Huvitz's outstanding optical technology and innovative image software. Shows extensive information, such as 3D structure of Retina, Macula's thickness and separation in a vivid image.

High Resolution Image - min. 60 lines/mm of central Fundus

Creates 3um OCT Digital Resolution medical images, allows more precise Retina observation and useful follow-up examinations.

Accurate and Stable Image Averaging

It is very important to obtain high-quality images that are accurate and stable in all OCTs. However, it's not easy to capture these due to patient eye movement over the period of the test. The HOCT detects fast eye movements with image processing algorithms of fast Scan Speed* and Smart Viewing Technology(SVT)** and scans up to 68,000 points per second, and calibrates to create a high-quality optical image.

HOCT can acquire high-quality images without any repetitive operation for first time users.

* 68,000 A-scan / sec., Less than 1.4 sec. in 6x6 mm² 3D shooting

** Smart Viewing Technology : Huvitz's Speckle-Noise-Reduction System & Pre-Acquiring Algorithm to acquire high-quality images


Vividly Visualized Retinal Layers

Visualizing with precise B scans and smooth 3D Images at faster scan speed makes it easier to observe pathological shapes and status in stratified Retinal Layers.

It is also useful to further elucidate the pathological rheobase of Macula and Optic Disc, including factors that impair Photoreceptor Function, Retinal & Choroidal Vasculature(vascular system) in a slice image for Retinal Layer consists of 7 pieces.

Brightness Level Adjustment

Precisely identify lesions by minutely adjusting image's brightness and contrast. In this way, specific parts of lesions can be highlighted which help users to easily see details.

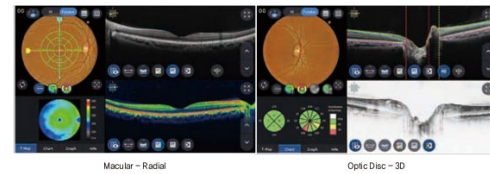



One for All System

3D OCT Angiography, Fundus Camera, Built-in PC :
 Combined-One is more accurate and useful

By combining OCT Angiography, Full Color Fundus Camera, and PC, it can generate high resolution images providing multi-purpose functions for diagnosis. It saves both time and space by performing frontal view(Enface) of eye diseases, Tomography, cross-compare and diagnosis in a single run.

Upgraded Combined-One

Provides maximum psychological stability to the patient without re-shooting and reduces stress during shooting*. Easily checking lesion's position by Fundus Image, it precisely guides the location of the OCT Scan Image.

*Motion detection technology/Smart Scan Technology (SST) is applied to achieve perfect images without re-shooting even though there're flicker or movement (see Smart Scan page).

Compact Design - it can be installed in a small space

Thanks to HOCT's space-saving design, it is perfect for hospitals and research areas with a variety of diagnosis devices & treatment equipment. HOCT can maximize the convenience for users as well as patients by saving time & space.

Web Browsing System to view data anytime, anywhere

Patient's test data can be analyzed anywhere on the Internet. You can check and analyze all data of HOCT through Web Browser such as Internet Explorer, Safari, Chrome without installing special software separately.



User Friendly

Auto Tracking & Auto Shooting :
 makes it easy to use and obtain reliable data


HOCT is smart.- Obtains reliable data with minimum deviation of image quality according to user's measurement proficiency.

Fast and Stable Full Auto Mode

Simply press the button once to capture the image quickly and easily without any errors with Auto Tracking, Optimize, Auto Shooting at the correct position. Depending on the application, select Semi Auto to obtain more detailed images.

Semi Auto Mode for more precise images


You can obtain a more precise image by shooting Semi Auto Mode turning one's gaze to the side for patients with eye diseases such as cataract, strabismus, or optic disk and peripheral measurements. Semi Auto Mode can also be applied to eyes with weak signals. XY alignment, focus is automatically adjusted, and manual operation during auto adjustment is also possible. Focusing and Firing functions can be judged and involved by users so that users can obtain images in an intuitive way.

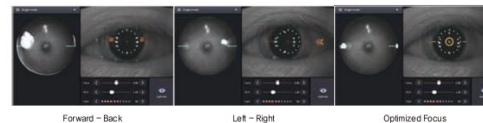

Smart Scan

Start and finish instantly through only one-click :
Speedy process reduces errors in forward looking of patients

It provides convenience & accuracy by offering easy & various scanning functions with Macula, Optic Disc, and Anterior.

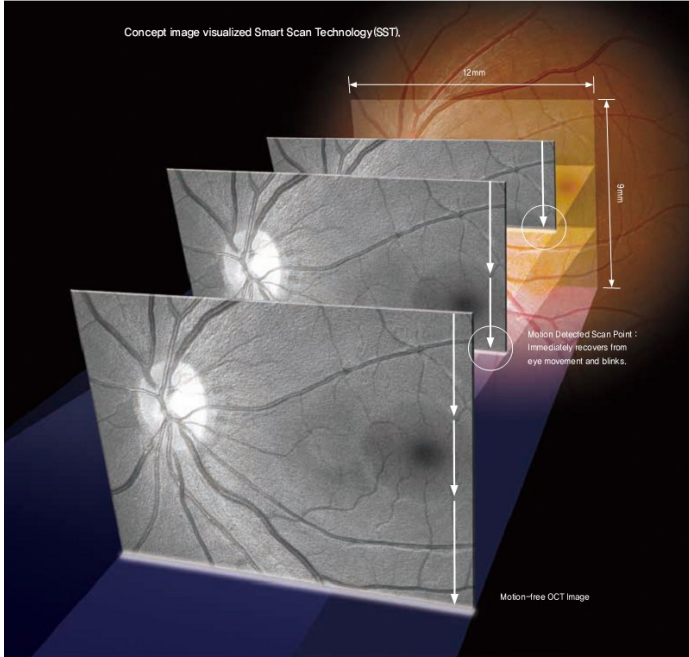

Wide Area Scan(12mm x 9mm) for efficient diagnosis


A quick scan covers Macula and Optic Disk areas extensively. By scanning around Optic Disc or Macula for patient's pathological status, you can check the Thickness Maps between RNFL(Retinal Nerve Fiber Layer), GCL(Ganglion Cell Layer) and RPE Layers.

Smart Scan Technology with Motion Detection

Image analyzer with Huvitz's unique Smart Scan Technology(SST) obtains a complete and perfect C-scan image by detecting any motion of eye flicker or movement that would prevent disappearance of scan line and image collection during measurement.

Providing Various and Useful Scan Patterns

12 different patterns make it available to choose and apply the optimized pattern to the main symptoms or the area of retinal disease without repetitive work or time-wasting.

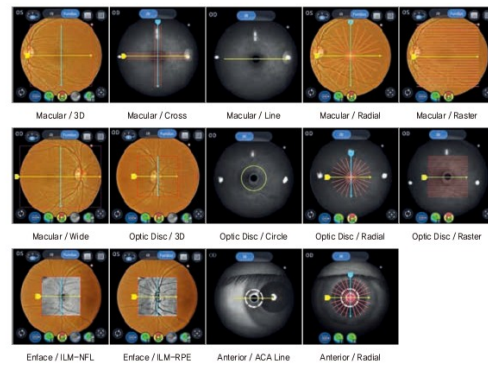

Accurate Analysis

Accurate segmentation and measurement :
Analyze pathology status from various perspectives

A complete analysis helps you observe symptoms, illnesses and progress of each patient at a glance. Key indicator values compared to Normative Data are displayed in table and chart format.

Progression to track pathological changes

OCT scan and fundus image of a patient can be compared at a glance to sequential measurement results from baseline to present.
Progression from past to present helps analyze disease progression and treatment process.
Thickness, Enface, and ETDRS can be superimposed on the IR or Fundus at each measurement point so that the change in thickness of nerve fiber can be confirmed according to the transition. It also provides a trace graph so you can study at a glance.

Compare before and after patient's symptom

You can compare and analyze the baseline data of a patient with the current data.

3D modeling in high speed and wide area


High-speed, wide area(12mm x 9mm) 3D images help you quickly and comprehensively understand the condition of the Retina. Also, layer thickness maps can be used from ILM to RPE, respectively and Morphological changes on the measured surface of the layers can also be visually confirmed.

OU to cross-analyze function of binocular

Provides comparative analysis for Macular Thickness, RNFL Thickness, ONH(Optic Nerve Head) of binocular.

Summary : Monocular-Scan and OCT / Fundus image


Provides a summary analysis of Macula retina, RNFL, ONH at a glance. Helps identify whether follow-up examinations are needed or not. Easy to explain the results to the patient after diagnosis.

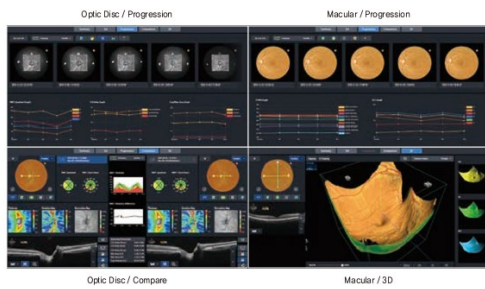

Detailed Report

From quick summary to simple comparison and complex evaluation : Complete a perfect report

Provides patient's pathological structure and relevant & important data in easy-to-read format and also can print out the report on analysis screen.
Analysis results can be viewed via Web Browser and printed out with different types of report.

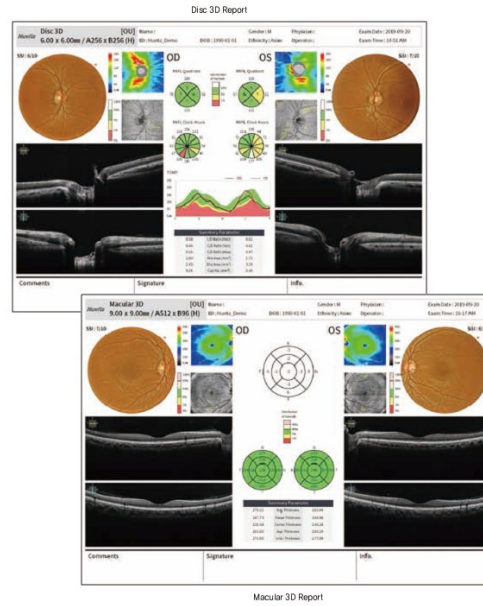



Anterior Measurement

One Single System :

Start and finish in one place, making patient more comfortable

Anterior Segment Module allows measurement and analysis of cornea thickness, angle and 3D image. It helps users work more efficiently by acquiring both anterior and posterior in one place.

9mm (16mm)* Wide Chamber View

Measurement of ACA(Anterior Chamber Angle) between cornea and Iris allows diagnosis and management of angle-closure glaucoma patients.

*9mm & 16mm of Anterior Segment Lenses are optional.

9mm High Resolution Cornea Thickness Measurement

9mm High Resolution Cornea Scan provides an objective view of the structure of the eyeball. It displays a cross-sectional image of the measured corneal thickness.

Corneal Thickness Map

Corneal's irregularity, Thinnest point can be identified with a corneal thickness map to visualize the patient's corneal thickness at a glance.



ACA Measurement

Corneal Thickness Measurement



Full Color Fundus Image

Insight of Posterior Segment of Eye : for Comprehensive diagnosis

Color Retinal Images optimized with high-resolution and contrast are very useful in analysis and clinical diagnosis.

Best images are provided by Low intensity of flash, fast capture speed, quiet operation, small pupil mode and automatic flicker detection.

High resolution and performance 12 Megapixel Camera

High performance camera with Motion Artifact Suppression Technique provides high resolution images and also its low intensity of flash, fast & quiet operation maximize measurement quality.

Auto-Detection Of Pupil Size and Auto Flash Level Function


It accurately measures the pupil size and automatically adjusts the intensity of light according to pupil size. Even patients with small pupil sizes can be easily measured without switching mode. Small Pupil Mode can be an option to adjust more intensive light for the small pupil size.

Fixation Target for flexible configuration

Fixation target can be set on the display for fine adjustment of a specific part of the eyeball.

Panorama function for wide range of peripherals

Multiple built-in capture color fundus images at different positions and automatically stitch them to optimized total overview. By providing high-resolution images with minimal distortion, you can immediately see key information for a comprehensive assessment of patient's eye.

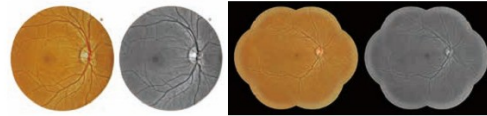

Fundus Image

Panoramic Image/Non-Mydiatic Composite Retinal Image

Innovative Angiography

Auto-Analysis of Retina & Microvessel of Choroid :
Customized Treatment per Patient with Details

Huvitz's own optical technologies, Real Time Tracking / Noise Cancelling / Motion Correction are inter-operating, and automatically analyzing & visualizing Retina, Microvessel of Choroid.

The Innovation of Optical Technologies :
Faster & More Convenient Analysis is Now Able

HOCT provides reliable results by minimizing measurement errors with combination of Huvitz's own optical technologies which are called as TAT(Triple Angiography Technology).

Three Innovative Optical Technologies – TAT(Triple Angiography Technology)

Real Time Tracking


With A-scan 68,000/per second, Real Time Tracking technology creates connecting vessel images by minimizing Motion Artifact caused by eye blinking, movement of eyeballs.

Noise Cancelling




Visual Processing Algorithm refines details of images in an instant. Therefore, it visualizes vessel conditions of Retina's layers as high quality images.


Motion Correction

Motion Correction Technology helps to make the Angiography images without distortion by line correction of vessels' dislocation.

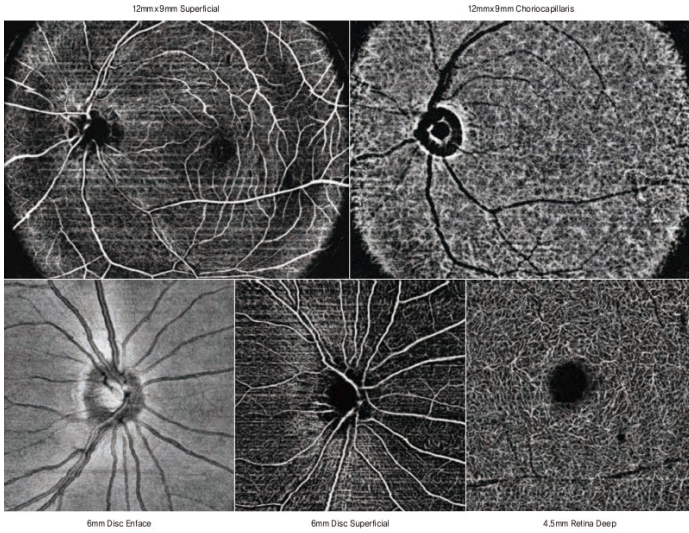

Real Time Tracking Technology Algorithm

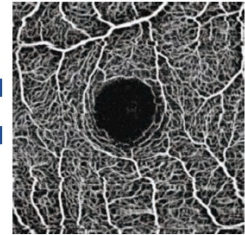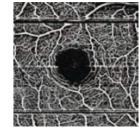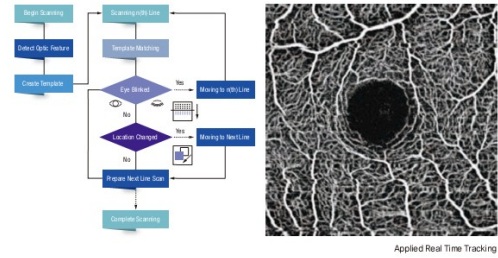

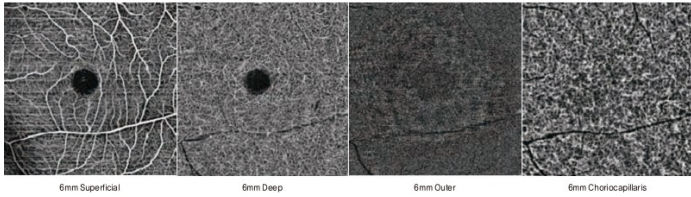

Innovative Angiography

68,000 scan/sec by One Button :
High-Resolution Images with Quantified Index


By One-Button, accurate details are provided with high-resolution images for vessels of retina & choroid and data of FAZ, flows, density.

Retina Layer Auto-Analysis

The function provides high-resolution images for vessels of retina & choroid with quantified index. It can be utilized with early diagnosis & treatment progression for macular degeneration, diabetic retinopathy, glaucoma, hypertensive retinopathy and retinal vein occlusion. Also, users can check abnormal vasculature in Custom View. Since the analysis shows index for FAZ, Flows, Density, it's easy to establish treatment plan.

Detail Display for Accurate Index and Evaluation

In Detail Mode, users can specifically observe vascular network per layer. Using analyzing tool, details of FAZ can be acquired conveniently.

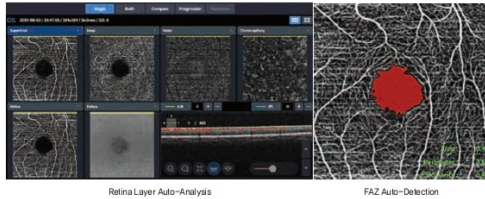

Retina Layer Auto-Analysis


FAZ Auto-Detection

Progression

The Progression mode helps users to follow up the pathology of a disease.

Binocular Comparison(OU)

In the Comparison Mode, users can check vascular network by layer in detail. By indicating layers in different colors, it's easy to check and understand the pathology of a disease. In case of comparison for diabetic retinopathy, the mode helps to track the pathology and establish treatment plan.

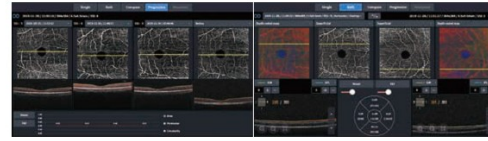


Progression

Binocular Comparison

A Variety of Scan Sizes : 3mm / 4.5mm / 6mm / 9mm

HOCT-Angiography supports a variety of scan sizes, users can choose and observe per needs & cases.

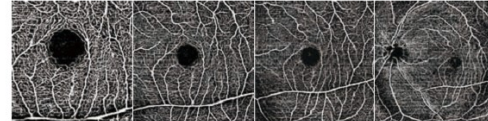

Various Scan Sizes : 3mm / 4.5mm / 6mm / 9mm

Angio Panorama

In case of checking Angiography image with large size, it is convenient to utilize the Angio Panorama function.

*3mm(Max. 12mmx9mm), 4.5mm(Max. 13.5mmx9mm) & manual mode are available




Panorama Mode / Panoramic Image

1TB Storage


By internal 1TB Storage of HDD, users can manage data with enough space.

Clinic Exams


High-quality, high-resolution OCT and color fundus images from HOCT are extremely useful for analysis and clinical diagnosis as the pathologic structure and status of each layer are accurately observed and recorded.

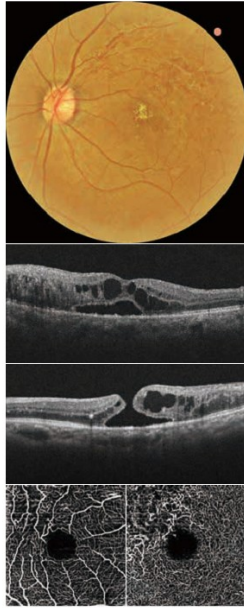


Macular Hole (MH), RVO (Retinal Vein Occlusion)
 A macular hole is a retinal break commonly involving the fovea. Severe stage of RVO progressed to a macular hole.

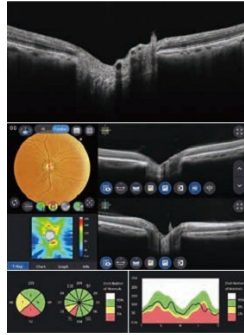


Glaucoma
 Glaucoma is a disease that damages your eye's optic nerve. The same symptoms are found at Thickness map, Fundus, TSNIT chart, Clock chart.

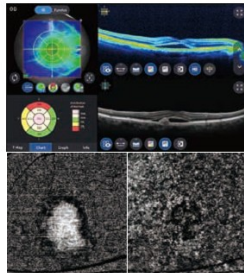


Macular Degeneration(MD)
 Age-related macular degeneration is a disease that blurs the sharp, central vision you need for straight-ahead activities.

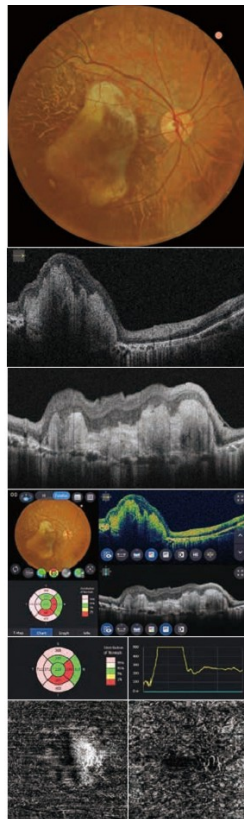

Choroidal Neovascularization (CNV)
 Choroidal neovascularization refers to condition that new blood vessels created in the choroid layer. In this case, because of huge choroidal neovascularization, distort the normal alignment of the overlying retina. It can be seen that the overlying photoreceptor layer is not clearly delineated from the other retinal layers, and the temporal retina also has the geographic atrophy of the photoreceptor layers and retinal pigment epithelial layers. This results is also visible on the Fundus photograph and Outer retina of angiography.

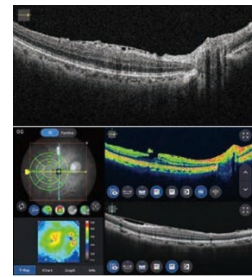

Epiretinal Membrane (ERM)
 Epiretinal membrane is a disease of the eye in response to changes in the vitreous humor or more rarely, diabetes.

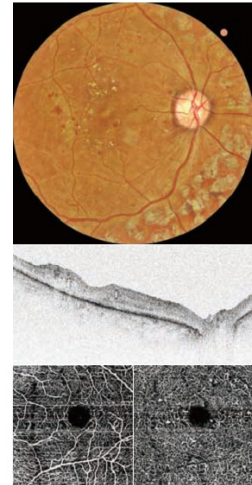

Diabetic Retinopathy (DR)
 Diabetic retinopathy is when high blood sugar levels cause damage to blood vessels in the retina. These blood vessels can swell and leak. Or they can close, stopping blood from passing through. Sometimes abnormal new blood vessels grow on the retina.