

VetiGen P50

Veterinary Portable X-ray Unit

Innovate Ahead, Lead Simplicity

companion animal



Q | Search Mindray Animal Medical



LinkedIn Facebook Instagram


www.mindrayanimal.com
service@mindrayanimal.com

©2026 Shenzhen Mindray Animal Medical Technology Co., LTD. All rights reserved.
P/N:ENG-VetiGen P50-210285X4P-20260408

Easier Operation, Higher Efficiency

iShotHelper

Professional positioning guidance, friendly to new vets, greatly reduce the operation threshold.

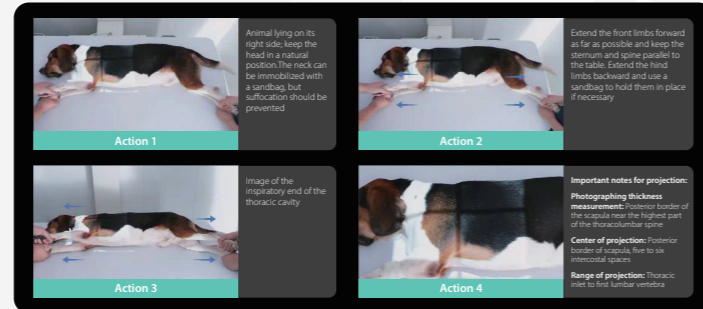

Precise Tilt & Angle Display

The machine tilts freely, with real-time angle shown on screen for accurate shooting setup.

Foot Switch & Hand Switch

- Standard equipped with hand switch & foot switch dual control.
- Enabling single-person to complete the shooting and greatly improving shooting efficiency.



Touchscreen

- 10"capacitive touchscreen for simple operation.
- Designed with inclination for easy observation without glare interference.

Professional Design, Safe Protection


APR

- Built-in cat/dog APR, with customizable settings for exotic pets.
- Automatically sets mA/ms and mAs based on body thickness.




Multi-color Status Indicators

The exposure indicator light displays multiple states, clearly distinguishing the three states of start-up, preparation and exposure, making the operation process safer.

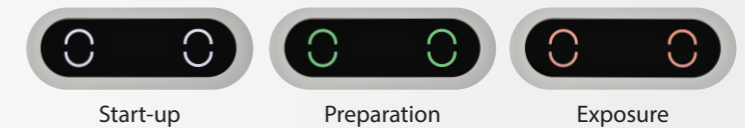

Silent Operation

Silent operation for animal comfort and reduced stress, improving the shooting experience.

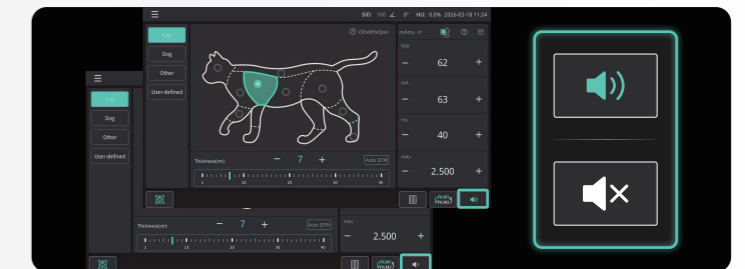

Technical Specifications

| Technical Specifications | |
|----------------------------|--|
| Construction | Monoblock unit with touch screen High-frequency technology |
| Max Power Output (kW) | 5 kW |
| Tube Voltage (kV) | 40~125 kV, 1 kV per step |
| Tube Current (mA) | 12.5~100 mA |
| Current time product (mAs) | 0.16~160 mAs |
| Loading time (ms) | 1.6 ms~8000 ms |
| Working Frequency | 400 kHz |
| Focus spot (mm) | 0.6 mm / 1.8 mm |
| Total Filtration | 2.9mm Al collimator included |
| Size (L/W/H) | 395 x 285 x 245 mm |
| Weight (kg) | 15.8±0.2 kg |
| Foot switch cable (m) | 3 m |
| AUTO-BTM | Option |